

Chat Editor: How to Import Pages from a Different Vocabulary File into Your Student's Vocabulary File

Introduction

In this tutorial, you will learn how to import different pages from a variety of vocabulary files within TouchChat directly into your student's custom vocabulary file. This way you can save time from re-inventing the wheel and use a variety of already made templates that are easily accessible for custom tailoring to your students learning needs.

Learning Objectives

Completion of this tutorial will give you experience with how to:

- Import another page from a different page set.
- Locate the newly imported page within your vocabulary file.
- Link a page to an existing button.

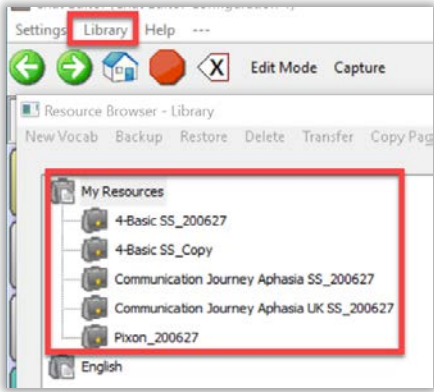
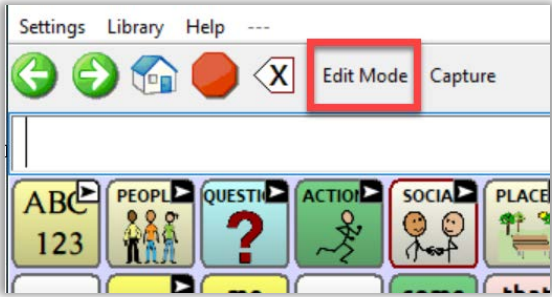
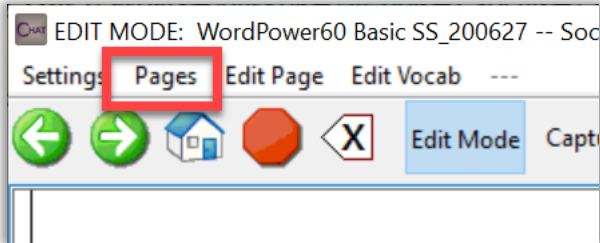
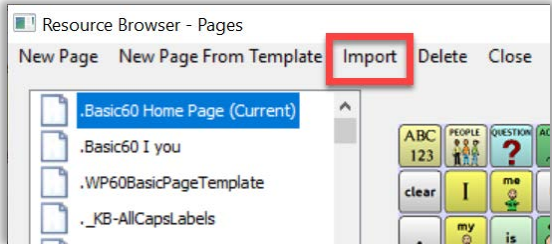
This tutorial assumes that you have

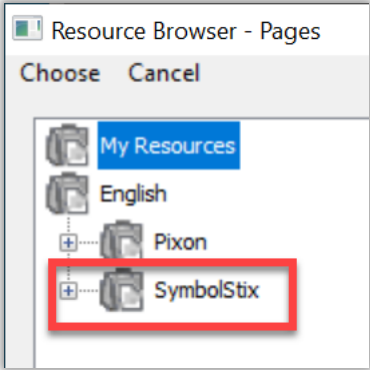
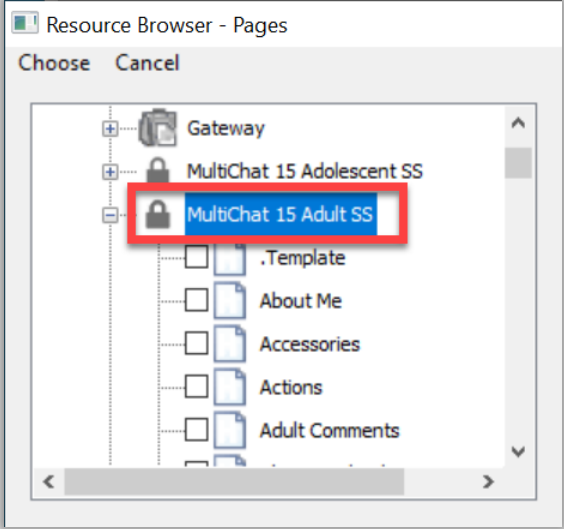
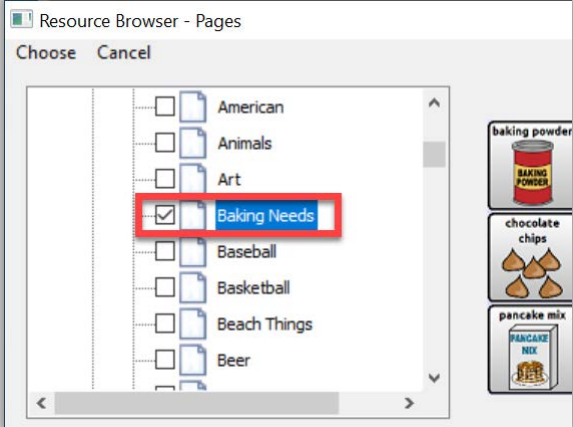
- An iPad with TouchChat installed, a NovaChat, or Lamp Words for Life.
- The most recent version of Chat Editor installed on your PC.

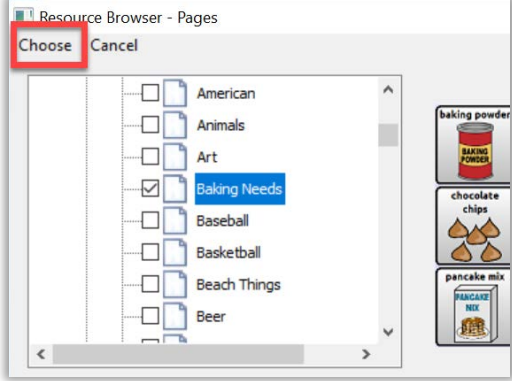
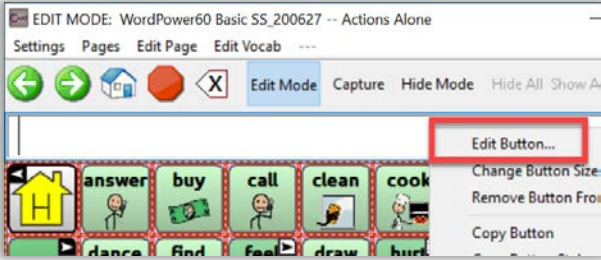
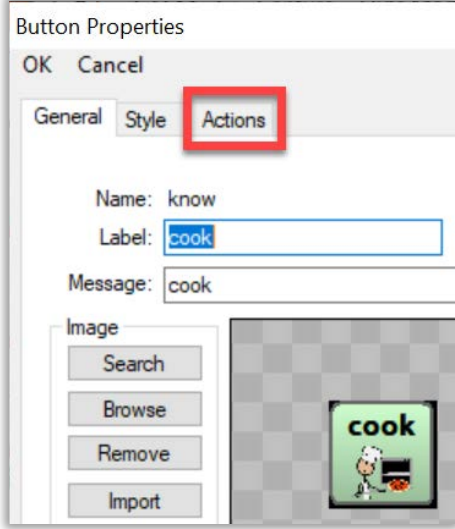
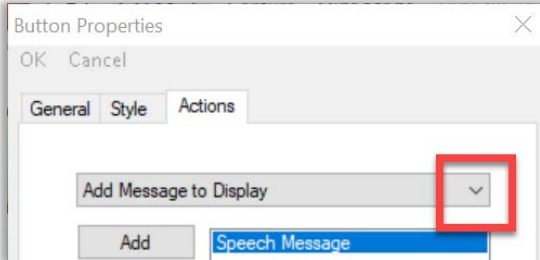
Case Study

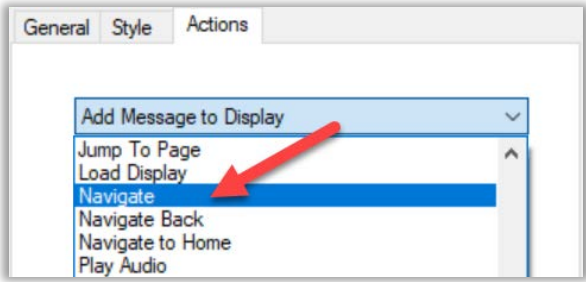
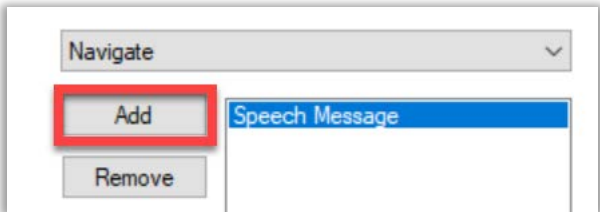
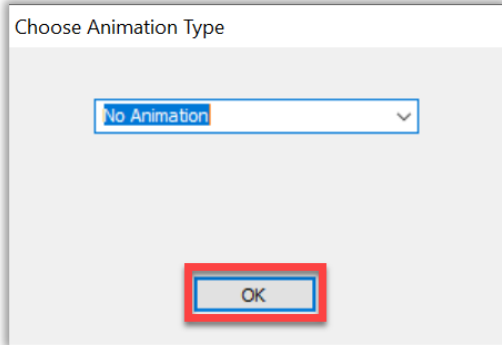
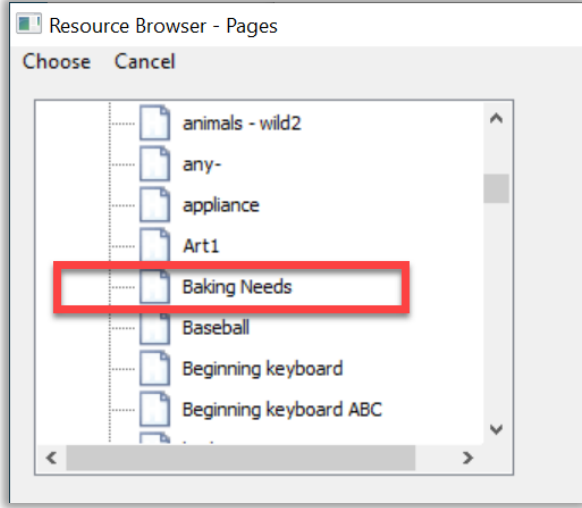
Cassie is in grade 2 and she recently received a SET-BC iPad with TouchChat. Her team is excited to customize the WP40 pageset for Cassie. After initial training with a SET-BC consultant on how to use TouchChat and its many different vocabulary files and pre-made templates, her team likes the look of some of the activity pages on Multi-Chat 15. Rebecca, Cassie's EA, knows that Cassie is really motivated with baking and saw a wonderful template in Multi-Chat 15 she hopes to embed in Cassie's original WordPower 40 pageset. Now, she will follow this tutorial to import pages from a different pageset.

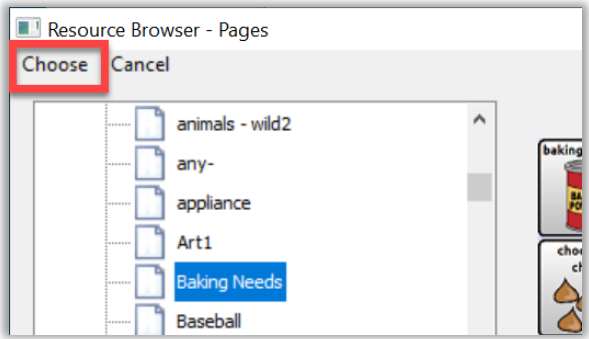
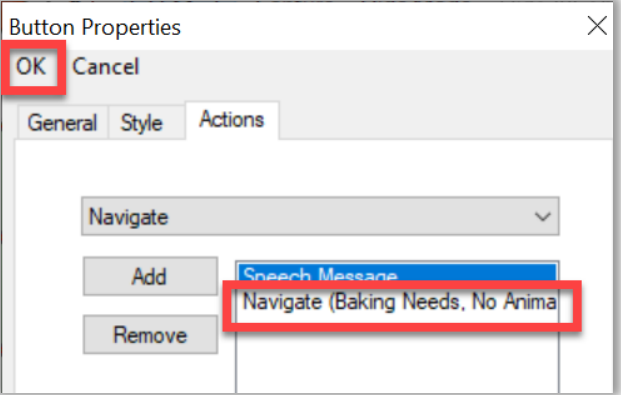
How to Import Pages from a Different Vocabulary File into Your Student's Vocabulary File

<p>1</p>	<p>Launch Chat Editor.</p> <p>To import an existing page to your student's vocabulary set, first open the students' vocabulary file.</p> <p>Click on Library, and you should be able to locate the file under the My Resources folder.</p>	
<p>2</p>	<p>Select Edit Mode.</p> <p>Then, determine the button which the imported page will be linked to.</p> <p><i>If you would like to create a different button, refer to Module 6 on how to customize buttons.</i></p>	
<p>3</p>	<p>To find vocabulary templates to be imported, select Pages.</p>	
<p>4</p>	<p>Click on Import.</p>	

<p>5</p>	<p>You will be able to select a vocabulary page from the My Resources or the English folders.</p> <p>The My Resources folder is where you may have backed up students' vocabulary files using Dropbox.</p> <p>To locate preset TouchChat and Nova Chat vocabulary files, click on English, and then SymBolStix.</p>	
<p>6</p>	<p>For example, under SymbolStix, click on MultiChat 15 Adult SS to expand the folder.</p>	
<p>7</p>	<p>Click on the checkbox beside the vocabulary page you wish to select.</p> <p>For example, click on the checkbox beside Baking Needs.</p>	

<p>8</p>	<p>Click Choose on the toolbar.</p> <p>Note that more than one page can be selected to import.</p> <p>The chosen pages are now imported to the student's vocabulary file.</p>	
<p>9</p>	<p>To link the imported page to the button determined in Step 2, right click on the button.</p> <p>Then, select Edit Button.</p>	
<p>10</p>	<p>Select the Actions tab.</p>	
<p>11</p>	<p>Click on the downward arrow next to "Add Message to Display" to expand the dropdown menu for different options.</p>	

12	Select Navigate .	
13	Then, click on Add .	
14	Select No Animation and click OK .	
15	All the pages in this vocabulary file will be loaded. Locate your imported file from the list.	

<p>16</p>	<p>Click on Choose on the toolbar.</p>	
<p>17</p>	<p>Select the imported page and then click OK.</p>	
<p>18</p>	<p>You have now linked the imported vocabulary page to the button.</p>	